

SUPAPORE TP

AMAZON

PLEATED MEMBRANE FILTER CARTRIDGE

SUPAPORE TP Microfiltration cartridges have been designed specifically for use in a wide range of aggressive gas and liquid processing applications. Its naturally hydrophobic PTFE membrane makes this cartridge an ideal choice for systems requiring excellent compatibility, reliability and a consistently high quality filtrate.

Product Features

SUPAPORE TP cartridges feature:

- A naturally hydrophobic PTFE membrane rated at 0.2µm in liquids and 0.01µm in gases
- Polypropylene support
- Minimal fibre migration
- Wide chemical compatibility
- Validation by bacterial and aerosolised bacteriophage challenge

All cartridges are thermally bonded and adhesive free. Each unit is pre-flushed with ultra-pure water and integrity tested before final assembly*1.

- Can be steam sterilised
- Sanitisable
- Meet current US FDA CFR Title 21 requirements
- Fully integrity testable by Water Intrusion, Bubble Point or Diffusion Flow test*1
- All B grade modules are integrity tested during manufacture

Available in a wide range of grades & absolute rated pore sizes:

- 0.1, 0.2, 0.45 and 1.0µm
- General Grade
- Biological Grade

Suitable for most sanitisation regimes including steaming, autoclaving, hot-water flush and most common sanitising agents.

(Note: Steaming, autoclaving and hot water sanitisation must use the stainless steel support ring or the glass filled end-cap option)

Amazon Filters manufactures a comprehensive range of filter housings for **SUPAPORE TP** cartridges including the industrial **50** and **60 Series** and **the 70 Series** for hygienic applications. Please consult our Sales Office for further details on the full range.



SUPAPORE TP FEATURES AND BENEFITS

- Validated*2 removal ratings for consistent and reliable performance
- Fully Integrity Testable - Water Intrusion, Bubble Point and Diffusional Flow values provided for all biological grade pore sizes*1
- Feature naturally hydrophobic PTFE membrane and Polypropylene hardware
- Thermal bonding process eliminates adhesives and ensures minimal extractables
- Identification mark on every cartridge
- Full product guide available

*1 B grade only

*2 See product guides for further details stating test method and rated efficiencies

SUPAPORE TP

SUPAPORE TP

TECHNICAL DATA

Dimensions

Outside Diameter: 70mm
 Typical Surface Area: 0.73-0.80m² (Per 10")
 Length: See ordering guide

Sterilisation and Sanitisation*³

Steam or Autoclave: 136°C for 30 mins (100 cycles)
 Hot Water: 90°C for 30 mins (0.2 bar Δp max)

*³Applies to single open end cartridges only. For all steaming and hot water applications, the stainless steel end cap insert or Glass Filled end cap option must be used.

Maximum Operating Conditions

Temperature: 80°C

Recommended Maximum Differential Pressure:

Forward Flow: 5.0 bar @ 20°C

Reverse Flow: 2.1 bar @ 20°C

Product validation guide available on request. All SUPAPORE TP cartridges are manufactured under strict control with batch number identification, giving full traceability on all components.



Materials of Construction

Filter Media: PTFE
 Media Support: Polypropylene (PP)
 Cage, Core & End Caps: Polypropylene (PP)

Integrity Test Specification

B Grade

Micron Rating	Bubble Point In 60/40 IPA/Water	Diffusion Flow (cc/min/10") In 60/40 IPA/Water
0.10µm	≥ 1620mbar (23.5 psi)	≤ 10 @ 1000mbar
0.22µm	≥ 1080mbar (15.7 psi)	≤ 10 @ 800mbar
0.45µm	≥ 630mbar (9.1 psi)	≤ 10 @ 400mbar

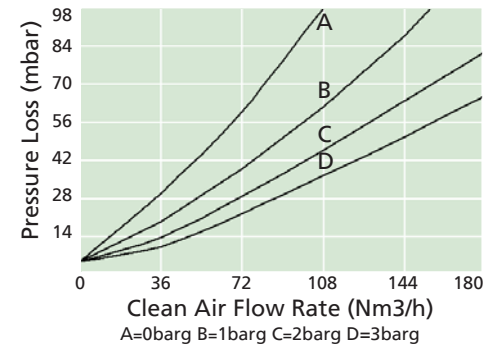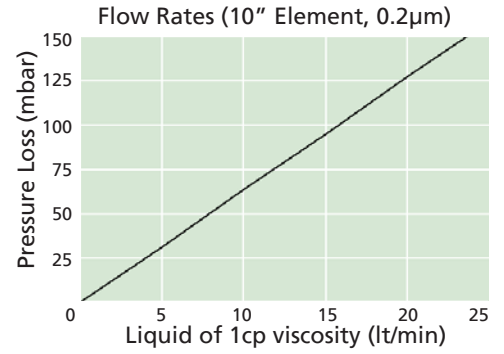

INDUSTRIES AND APPLICATIONS

Food and Beverage
Pharmaceutical
Fine Chemicals
Electronics

- CO₂ lines, Sterile vent for storage and transfer tanks
- Solvent filtration, polishing
- Polymers, Aggressive acids, Bases, Solvents, Ozonated DI Water
- Photoresists, Acids, Etch Solutions

AMAZON FILTERS LTD.

ALBANY PARK ESTATE, FRIMLEY ROAD,
 CAMBERLEY, SURREY, GU16 7PG
 ENGLAND

Tel: +44 (0) 1276 670 600

Fax: +44 (0) 1276 670 101

Email: sales@amazonfilters.co.uk

Web: www.amazonfilters.co.uk



Certificate No. FM39584

Distributed by: